

HD Digital Video Colposcope



*Don't take a chance
with Diagnosis
Perfect Picture
For Perfect Diagnosis*

MODEL : KM-20K HD

Technical Specifications

- 1/6" Microlens SONY HD CCD
- Imaging System from Original SONY
- Resolution >950 CCD Lines
- Digital Matrix Processor Technology
- Effective 20 Mega Pixels
- Multiple modes of Electronic Green Filter
- Black Morphological Filtering
- Automatic Screening, Diagnosis & Reporting
- Optical Zooming: 1-52x
- Corner to Corner Uniform Brightness
- Output : RS-232, S-video, composite and BNC
- Remote operated
- No need of dongle and capture card

KINESIO MEDICARE

- Designed according to the guidelines of IFC.
- Integrated Management of View Examination, Image Capture & Freeze, Recording.
- Can be integrated to LAN and HIS
- Reference Picture library to help user.
- Statistics mode with automatic flow chart facility is provided.
- Printing of Multi Format photo report with Hammond graph Report
- Powerful Comparison to compare.
- Colposcopy REID Evaluation.
- Report can be e-mailed as well.

KINESIO MEDICARE



ISO 9001:2015

Admin Add: 1-54/14 Off No. 201, Sai Durga Residency, Ramakrishna Nagar, Hyderabad - 500049

Contact Us: 9618133196, 9100028443

Ultra HD Digital Video Colposcope



*Don't take a chance
with Diagnosis
Perfect Picture
For Perfect Diagnosis*

MODEL :KM-40K Ultra HD

Technical Specifications

- Full HD Technology with CMOS Sensor
- 1080p Resolution Sony Sensor (1920x1080)
- Aspect Ratio 16:9
- 21.2 Mega Pixels
- No need of Dongle & Capture Card
- Multiple Grades of green Filter
- Recommended Lamp life 60000 hours
- Highest quality-1-52x Magnification
- Acetic acid Test Timer
- Latest Digital Matrix Processor
- Fine Advance & fast Auto/Manual focus - Fine Technology
- Black Morphological Filter
- Video output : HDMI, Composite and many more
- Remote operated

KINESIO MEDICARE Image Processing Software

- Designed according to the guidelines of IFC.
- Integrated Management of View Examination, Image Capture & Freeze, Recording.
- Can be integrated to LAN and HIS
- Reference Picture library to help user.
- Statistics mode with automatic flow chart facility is provided.
- Printing of Multi Format photo report with Hammond graph Report
- Powerful Comparison to compare.
- Colposcopy REID Evaluation.
- Report can be e-mailed as well.

KINESIO MEDICARE



ISO 9001:2015

Admin Add: 1-54/14 Off No. 201, Sai Durga Residency, Ramakrishna Nagar, Hyderabad - 500049

Contact Us: 9618133196, 9100028443



Hysteroscope

- Ultra HD 4mm and 2.9mm 30° scope.
- HD vision with wide angle.
- Fully Autoclavable.
- Made in Germany.


Operative Sheath

- Quick lock locking system.
- Medical grade German steel.
- Size : 2.9mm and 4mm
- Instrument channel of 5fr.



TCRE Set

- Bi-polar working element with inner and outer sheath.
- Continuous flow outer sheath.
- Compatible with all German scopes.
- Size : 26fr for 4mm scope and 18.5fr for 2.9mm scope.

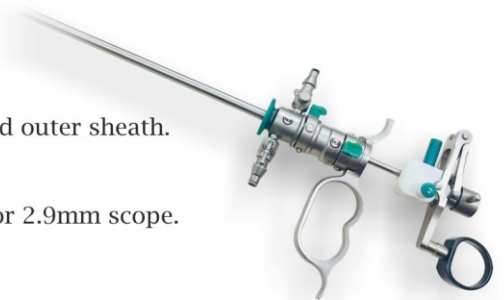

Hysteroscopy Pump

- Automatic German sensor operated.
- Maximum flow rate 1 ltr/min.
- Tissue sensing technology.
- Maximum pressure upto 400mm Hg.



Hysteroscopy Instruments



Grasper



Scissor



Biopsy

KINESIO MEDICARE



ISO 9001:2015

Admin Add: 1-54/14 Off No. 201, Sai Durga Residency, Ramakrishna Nagar, Hyderabad - 500049

Contact Us: 9618133196, 9100028443

Authorised Distributor
TFS-UK



BLUETOOTH MODEL (SVDS-BT)



Contact Us

+91 9618133196 | +91 9100028443



Our Mail

kinesiomedicare@gmail.com



Our Address



1-54/14 Off No. 201, Hyderabad

VACUUM DELIVERY SYSTEM



MODEL
DIGITAL LED
MODEL WITH
ALARMING
(SVDS-02)

ANALOG MODEL

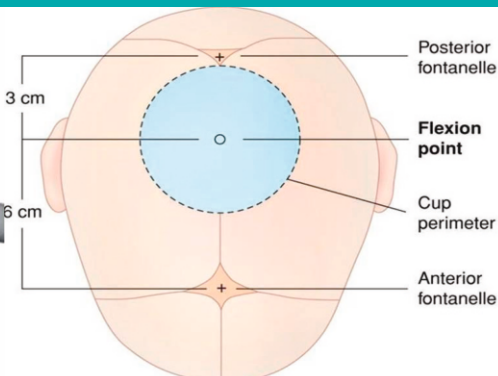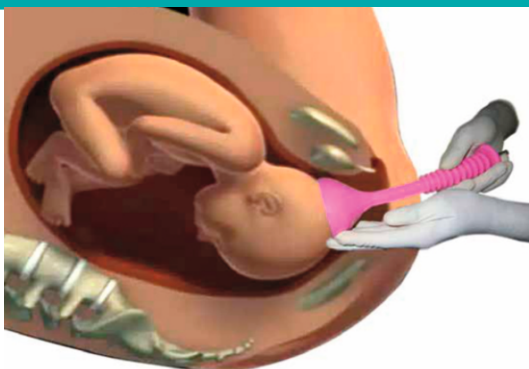

Kinesio SOFT SILICONE CUP :

- Medical-grade silicone rubber ensures absolute bio-compatibility and high durability
- Available sizes: 50 mm, 55 mm, 60 mm Suitable for occipito - anterior (OA) positions and simple outlet presentations
- Sterilization : Autoclavable
- Also we have birds metal cup and malmstorm design cup
- Sizes : 40 mm , 50 mm , 60 mm

FLEXION POINT LOCATION :

- The flexion point is situated on the sagittal suture 3 cm in front of the posterior fontanelle
- The flexion point may be located during vaginal examination by identifying the posterior fontanelle and then moving the finger anteriorly at a distance of approximately 3 cm along the sagittal suture
- The tip of the finger will mark the flexion point.
- If the cup is placed exactly on the flexion point, it assists the delivery of the smallest diameters of the fetal head.

Kinesio VACUUM EXTRACTOR :

- Kinesio Vacuum Delivery Machines are available in Two models. One is the Analog model (SVDS-01) and the advanced model is the Digital model (SVDS-02) with stainless steel trolley.

Features of SIDRA vacuum extractor :

- The innovative design of micro-processor controlled Digital display (LED) with precise indication up to 600mmHg accuracy which may not be seen in the conventional Analog bourdon vacuum gauge.
- Vacuum indication can be pre-selected in the digital model either in mmHg or in Kg/cm.
- The vacuum pressure doesn't drop due to the non-return valve and the doctor can concentrate only on delivery.
- Unbreakable plastic jar for easy handling, cleaning and re-assembling. The jar will help Create the vacuum and a small amount of fluid collection. This equipment is electrically operated and less than (40 W) . It can be run with inverter also . Inside the machine, the pump is oil less, noise less and maintenance-free.
- This equipment along with the stainless steel trolley. Hanging provision offer to foot switch.